ImpactTB/BAA: Standard Operating Procedures for Data Analysis

Colorado State University Coding Team

2022-08-08

Contents

1	Overview 1.1 About the project							
2	Experimental metadata	7						
3	Initial mouse characteristics							
4	Mouse Weights	11						
5	Colony forming units to determine bacterial counts	13						
	5.1 Data description	. 13						
	5.2 Read in data	. 14						
	5.3 Example one	. 16						
	5.4 Exploratory analysis and quality checks	. 16						
	5.5 Exploratory analysis	. 16						
	5.6 Identify a good dilution for each sample $\dots \dots \dots$. 17						
	5.7 Calculate CFUs from best dilution/Estimate bacterial load for							
	each sample based on good dilution	. 17						
	5.8 Create initial report information for these data	. 18						
	5.9 Sample ANOVA							
	5.10 Save processed data to database	. 20						
	5.11 Example two	. 20						
6	Enzyme-linked immunosorbest assay (ELISA)	21						
	6.1 Read in data from excel file	. 22						
	6.2 #Join the OD table with the information table	. 26						
	6.3 ELISA data analysis optimization							
	6.4 Use the data analysis models in our data	. 31						
	6.5 Creating a function for fitting model							
7	Flow cytometry	39						

4 CONTENTS

Overview

1.1 About the project

The objective of the Immune Mechanisms of Protection against Mycobacterium tuberculosis (IMPAc-TB) program is to get a thorough understanding of the immune responses necessary to avoid initial infection with $Mycobacterium\ tuberculosis\ (Mtb)$, formation of latent infection, and progression to active TB illness. To achieve these goals, the National Institute of Allergy and Infectious Diseases awarded substantial funding and established multidisciplinary research teams that will analyze immune responses against Mtb in animal models (mice, guinea pigs, and non-human primates) and humans, as well as immune responses elicited by promising vaccine candidates. The contract awards establish and give up to seven years of assistance for IMPAc-TB Centers to explain the immune responses required for Mtb infection protection.

The seven centers that are part of the study are (in alphabetical order):

- 1. Colorado State University
- 2. Harvard T.H. Chan School of Public Health
- 3. Seattle Children Hospital
- 4. [more]

Colorado State University Team and role of each member:

- Dr. Marcela Henao-Tamayo: Principal Investigator
- Dr. Brendan Podell: Principal Investigator
- Dr. Andres Obregon-Henao: Research Scientist-III
- Dr. Taru S. Dutt: Research Scientist-I
- [more]

1.2 About this book

The aim of this book is to provide data protocols and data collection templates for key types of data that are collected over the course of this project. By using standard templates to record data, as well as starting from defined pipelines to process and analyze the data, we aim to standardize the collection and processing of data across this project.

Here, we have built a comprehensive guide to wet lab data collection, sample processing, and computational tool creation for robust and efficient data analysis and dissemination.

Experimental metadata

Metadata for an experiment:

- species
- start_date
- end_date
- experimental_groups

Initial mouse characteristics

At the start of each experiment with a mouse model, we record several measurements or characteristics of each mouse. We record these measurements along with an identifier for each mouse (for example, based on tags or ear notches), so that we can later link the initial characteristics of each mouse with later measurements on the same mouse.

The values that we initially record for each mouse include:

- group: An identifier for the experimental group to which the mouse is assigned (e.g., "Control", "Group 1")
- group_detail: A longer description of the mouse's treatment group (e.g., "Vaccinated with vaccine candidate A at 4 and 8 weeks")
- notch_id: The number that corresponds with the ear notch pattern of the mouse, which together with the mouse's experimental group provides a unique identifier for each mouse in the experiment
- mouse_number: The mouse's group number and notch ID combined; this provides a unique identifier for each mouse in the experiment
- cage_number: The number of the cage to which the mouse is first assigned. This may change over the course of the experiment, as mice might be removed from a cage due to fighting, etc. Any of these later changes of cage will be recorded [where]
- sex: Whether the mouse is male ("m") or female ("f")
- age:
- strain

We have created a spreadsheet template that can be used to record these data, which you can download by clicking here. Here is an example of how the first rows of this template might look once it's filled out:

This template should be used at the initial time when mice are brought into the experiment. The file format is an Excel file, so you can use it by saving it to your computer and then opening and recording data with Excel. Later code in this chapter will read in a file in this template format to provide basic summaries of the data. Later code will read in these files to record the data in a project-wide database, which will allow us to integrate it with other data collected over the course of the experiment.

[Rules for naming the file. Include experiment name / study ID?]

Mouse Weights

4.0.1 Overview

Extreme weight loss and loss of muscle mass, also known as cachexia, typically presents along side chronic inflammatory illnesses like Tuberculosis disease @ [baazim2022interplay]. We now recognize that cachexia is part of a systemic response to inflammation, and has been linked to upregulation of pro-inflammatory cytokines such as TNF, IL-6, and IFNg in humans @ [baazim2022interplay]. Additionally, studies support the role of key immune cell populations such as CD8+ T-cells that, when depleted, counteract muscle and fat deterioration @[baazim2019cd8], and suggest that CD8+ T-cells may metabolically reprogram adipose tissue.

In recognition of cachexia related illnesses and diseases, we tracked the progression of weight loss over the course of this study, as is done with many TB-mouse studies @ [smith2022host], @ [segueni2016controlled]. These data is also useful when correlating to CFU count as well as expression of cytokines and other biological markers @ [smith2022host]. Here, mice are weighed in grams weekly to monitor clinical status as TB patients frequently display weight loss as clinical symptom associated with disease progression.

The following contains information about how the data was collected, organized, and curated for analysis in RStudio.

4.0.2 Parameters

Weights are recorded in an excel worksheet.

Column titles are as follows: who_collected date_collected sex dob notch_id mouse_number weight unit cage_number group notes

Groups included are: bcg, saline, bcg+id93, saline+id93, saline+noMtb

The notes column contains information regarding clinical observations.

good reference: https://elifesciences.org/articles/74419#s4

```
library(readxl)
library(tidyverse)
```

```
## -- Attaching packages -----
                                          ----- tidyverse 1.3.1 --
                   v purrr
## v ggplot2 3.3.5
                            0.3.4
## v tibble 3.1.6
                   v dplyr
                            1.0.7
## v tidyr
           1.1.4
                   v stringr 1.4.0
## v readr
           2.1.1
                   v forcats 0.5.1
## -- Conflicts ----- tidyverse_conflicts() --
## x dplyr::filter() masks stats::filter()
## x dplyr::lag()
                  masks stats::lag()
```

4.0.3 Read in data

Data is stored in one excel sheet, each week is one sheet named as the date -> return vector for each sheet name

4.0.4 Can also use rio to read in the data, more streamlined

4.0.5 Clean data

```
dataset <- data$before_vaccination %>%
  select("sex", "mouse_number", "weight", "cage_number", "group")
# combining columns mouse_number and cage_number
dataset$mouse_id <- paste(dataset$mouse_number, "-", dataset$cage_number)</pre>
```

- 4.0.6 Body weight over time graph and statistics
- 4.0.7 Weight loss over time graph and statistics
- 4.0.8 Weight vs CFU
- 4.0.9 Weight vs ELISA results
- 4.0.10 Weight vs lesion burden

Colony forming units to determine bacterial counts

5.1 Data description

The data are collected in a spreadsheet with multiple sheets. The first sheet (named "[x]") is used to record some metadata for the experiment, while the following sheets are used to record CFUs counts from the plates used for samples from each organ, with one sheet per organ. For example, if you plated data from both the lung and spleen, there would be three sheets in the file: one with the metadata, one with the plate counts for the lung, and one with the plate counts for the spleen.

The metadata sheet is used to record information about the overall process of plating the data. Values from this sheet will be used in calculating the bacterial load in the original sample based on the CFU counts. This spreadsheet includes the following columns:

- organ: Include one row for each organ that was plated in the experiment. You should name the organ all in lowercase (e.g., "lung", "spleen"). You should use the same name to also name the sheet that records data for that organ for example, if you have rows in the metadata sheet for "lung" and "spleen", then you should have two other sheets in the file, one sheet named "lung" and one named "spleen", which you'll use to store the plate counts for each of those organs.
- prop_resuspended: In this column, give the proportion of that organ that was plated. For example, if you plated half the lung, then in the "lung" row of this spread sheet, you should put 0.5 in the prop_resuspended column.
- total_resuspended_uL: This column contains an original volume of tissue

- homogenate. For example, raw lung tissue is homogenized in 500 uL of PBS in a tube containing metal beads.
- og_aliquot_uL: 100 uL of th total_resuspended slurry would be considered an original aliquot and is used to perform serial dilutions.
- dilution_factor: Amount of the original stock solution that is present in the total solution, after dilution(s)
- plated_uL: Amount of suspension + diluent plated on section of solid agar

5.2 Read in data

```
library(readxl)
library(dplyr)
library(purrr)
library(tidyr)
library(stringr)
library(tidyverse)
library(gridExtra)
library(ggplot2)
library(ggpubr)
#Replace w/ path to CFU sheet
path <- c("DATA/Copy of baa_cfu_sheet.xlsx")</pre>
sheet_names <- excel_sheets(path)</pre>
sheet_names <- sheet_names[!sheet_names %in% c("metadata")]</pre>
merged_data <- list()</pre>
for(i in 1:length(sheet_names)){
  data <- read_excel(path, sheet = sheet_names[i]) %>%
    mutate(organ = paste0(sheet_names[i]))
  data <- data %>%
    #mutate(missing_col = NA) %>%
    mutate_if(is.double, as.numeric) %>%
    mutate_if(is.numeric, as.character) %>%
    pivot_longer(starts_with("dil_"), names_to = "dilution",
                  values to = "CFUs") %>%
    mutate(dilution = str_extract(dilution, "[0-9]+"),
           dilution = as.numeric(dilution))
```

```
merged_data[[i]] <- data
}
all_data <- bind_rows(merged_data, .id = "column_label") %>%
    select(-column_label)
head(merged_data)
## [[1]]
## # A tibble: 342 x 8
##
      count_date
                           who_plated who_counted groups mouse organ dilution CFUs
      <chr>
##
                                      <chr>
                                                  <chr> <chr> <chr>
                                                                         <dbl> <chr>
## 1 "\"February 21 2022~ BK
                                      BK
                                                                             O TNTC
                                                  group~ A
                                                               lung
## 2 "\"February 21 2022~ BK
                                      BK
                                                  group~ A
                                                               lung
                                                                             1 TNTC
## 3 "\"February 21 2022~ BK
                                      BK
                                                                             2 TNTC
                                                  group~ A
                                                               lung
## 4 "\"February 21 2022~ BK
                                      BK
                                                  group~ A
                                                               lung
                                                                             3 53
                                                                             4 9
## 5 "\"February 21 2022~ BK
                                      BK
                                                  group~ A
                                                               lung
   6 "\"February 21 2022~ BK
                                      BK
                                                               lung
                                                                             5 4
                                                  group~ A
## 7 "\"February 21 2022~ BK
                                      BK
                                                                             6 2
                                                  group~ A
                                                               lung
## 8 "\"February 21 2022~ BK
                                      BK
                                                                             7 1
                                                  group~ A
                                                                lung
                                                                             8 0
## 9 "\"February 21 2022~ BK
                                      BK
                                                  group~ A
                                                                lung
## 10 "\"February 21 2022~ BK
                                      BK
                                                  group~ B
                                                                lung
                                                                             O TNTC
## # ... with 332 more rows
##
## [[2]]
## # A tibble: 112 x 8
##
      count_date
                          who_plated who_counted groups mouse organ dilution CFUs
##
      <chr>
                          <chr>
                                     <chr>
                                                 <chr>
                                                         <chr> <chr>
                                                                         <dbl> <chr>
## 1 "\"April 25 2022\"" JR
                                                                             O TNTC
                                     JR
                                                 group_1 A
                                                                sple~
## 2 "\"April 25 2022\"" JR
                                     JR
                                                                             1 TNTC
                                                 group_1 A
                                                                sple~
## 3 "\"April 25 2022\"" JR
                                     JR
                                                 group 1 A
                                                                sple~
                                                                             2 53
## 4 "\"April 25 2022\"" JR
                                     JR
                                                 group_1 A
                                                                sple~
                                                                             3 9
##
   5 "\"April 25 2022\"" JR
                                     JR
                                                                sple~
                                                                             4 4
                                                 group_1 A
## 6 "\"April 25 2022\"" JR
                                     JR
                                                                            5 2
                                                 group_1 A
                                                                sple~
## 7 "\"April 25 2022\"" JR
                                     JR
                                                                sple~
                                                                             6 1
                                                 group_1 A
## 8 "\"April 25 2022\"" JR
                                     JR
                                                                             7 0
                                                 group_1 A
                                                                sple~
## 9 "\"April 25 2022\"" JR
                                     JR
                                                 group_1 B
                                                                sple~
                                                                             O TNTC
## 10 "\"April 25 2022\"" JR
                                     JR
                                                 group_1 B
                                                                sple~
                                                                             1 TNTC
## # ... with 102 more rows
head(all_data)
## # A tibble: 6 x 8
##
    count_date
                           who_plated who_counted groups mouse organ dilution CFUs
```

##		<chr></chr>			<chr></chr>	<chr></chr>	<chr></chr>	<chr>></chr>	<chr></chr>	<dbl></dbl>	<chr></chr>
##	1	"\"February	21	2022\~	BK	BK	group~	Α	lung	0	TNTC
##	2	"\"February	21	2022\~	BK	BK	group~	Α	lung	1	TNTC
##	3	"\"February	21	2022\~	BK	BK	group~	Α	lung	2	TNTC
##	4	"\"February	21	2022\~	BK	BK	group~	Α	lung	3	53
##	5	"\"February	21	2022\~	BK	BK	group~	Α	lung	4	9
##	6	"\"February	21	2022\~	BK	BK	group~	Α	lung	5	4

5.3 Example one

5.4 Exploratory analysis and quality checks

5.5 Exploratory analysis

Dimensions of input data:

##

<chr>

1 "\"February 21 2022\~ BK

Based on the input data, data were collected for the following organ or organs:

The following number of mice were included for each:

The following number of replicates were recorded at each count date for each experimental group:

The following number of dilutions and dilution level were recorded for each organ:

People who plated and collected the data. Date or dates of counting:

Based on the input data, the plates included in these data were counted by the following person or persons: Based on the input data, the plates included in these data were counted on the following date or dates:

```
all data %>%
  select(organ, who_plated, who_counted, count_date) %>%
  distinct()
## # A tibble: 3 x 4
     organ who_plated who_counted count_date
##
     <chr>
                                   <chr>>
            <chr>
                       <chr>
## 1 lung
            BK
                       BK
                                   "\"February 21 2022\""
## 2 lung
            BK
                       BK
                                   "\"April 18 2022\""
                                   "\"April 25 2022\""
## 3 spleen JR
                       JR
head(all_data)
## # A tibble: 6 x 8
##
     count_date
                           who_plated who_counted groups mouse organ dilution CFUs
```

<chr>

<chr>

BK

<chr> <chr> <chr>

lung

group~ A

<dbl> <chr>

O TNTC

## 2 "\"February 2	l 2022∖~ BK	BK	group~ A	lung	1 TNTC
## 3 "\"February 2	l 2022∖~ BK	BK	group~ A	lung	2 TNTC
## 4 "\"February 2	1 2022\~ BK	BK	group~ A	lung	3 53
## 5 "\"February 2	1 2022\~ BK	BK	group~ A	lung	4 9
## 6 "\"February 2	1 2022\~ BK	BK	group~ A	lung	5 4

Distribution of CFUs at each dilution:

Here's a plot that shows how many plates were too numerous to count at each dilution level:

Here is a plot that shows how the CFU counts were distributed by dilution level in the data:

5.6 Identify a good dilution for each sample

```
## # A tibble: 6 x 8
  count_date
                          who_plated who_counted groups mouse organ dilution CFUs
##
    <chr>
                          <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> 
                                                                       <dbl> <dbl>
## 1 "\"February 21 2022\~ BK
                                     BK
                                                 group~ A lung
                                                 group~ A lung
group~ C lung
## 2 "\"February 21 2022\~ BK
                                     BK
                                                                          4
                                                                                9
## 3 "\"February 21 2022\~ BK
                                    BK
                                                                          5
## 4 "\"February 21 2022\~ BK
                                                                          3 53
                                     BK
                                                 group~ D
                                                              lung
                                                                          2
                                                                                92
## 5 "\"February 21 2022\~ BK
                                     BK
                                                 group~ A
                                                              lung
## 6 "\"February 21 2022\~ BK
                                     BK
                                                                               7
                                                 group~ A
                                                              lung
```


5.7 Calculate CFUs from best dilution/Estimate bacterial load for each sample based on good dilution

```
# Calculating CFU/ml for every qualifying replicate between 10-75 CFUs. Column binding by organ meta <- read_excel(path, sheet = "metadata")
```

```
tidy_cfu_meta_joined <- inner_join(meta, tidy_cfu_data) %>%
  group_by(groups) %>%
  mutate(CFUs_per_ml = (CFUs * (dilution_factor^dilution) *
                          (total_resuspension_mL/volume_plated_ul) * 1000)) %>%
  select(organ, count_date, who_plated, who_counted, groups, mouse, dilution,
         CFUs, CFUs_per_ml) %>%
  ungroup()
## Joining, by = "organ"
head(tidy_cfu_meta_joined)
## # A tibble: 6 x 9
    organ count_date
                                who_plated who_counted groups mouse dilution CFUs
    <chr> <chr>
                                           <chr>
                                                       <chr> <chr>
                                                                       <dbl> <dbl>
## 1 lung "\"February 21 2022\~ BK
                                           BK
                                                                           3
                                                                                53
                                                       group~ A
## 2 lung "\"February 21 2022\~ BK
                                           BK
                                                                          4
                                                                                 9
                                                       group~ A
## 3 lung "\"February 21 2022\~ BK
                                                                         5
                                           BK
                                                       group~ C
## 4 lung "\"February 21 2022\~ BK
                                                                          3
                                           BK
                                                       group~ D
                                                                                53
## 5 lung "\"February 21 2022\~ BK
                                           BK
                                                                                92
                                                       group~ A
## 6 lung "\"February 21 2022\~ BK
                                           BK
                                                       group~ A
## # ... with 1 more variable: CFUs_per_ml <dbl>
```

5.8 Create initial report information for these data

```
tidy_lung_cfu_plot <- tidy_cfu_meta_joined %>%
  filter(organ == "lung") %>%
  mutate(group = fct_relevel(groups, "group_1", "group_2", "group_3", "group_4")) %>%
  ggplot(aes(x = groups, y = log10(CFUs_per_ml), fill = groups))+
  stat_boxplot( aes(x = groups, y = log10(CFUs_per_ml)),
      geom='errorbar', linetype=1, width=0.5)+
  geom_boxplot(aes(group = groups), fill = NA, show.legend = FALSE, color = "lightgrey
  geom_point(show.legend = FALSE)+
  labs(title = paste0("CFUs in early infected mouse lung"), x = "Group", y = "log10(CFTUs) color = "Group")+
  guides(shape = "none")+
  theme_minimal()+
  stat_compare_means(label = "p.signif", method = "t.test", ref.group = "group_1") +
  scale_y_continuous(expand = c(0, 0), limits = c(0, 8))
```

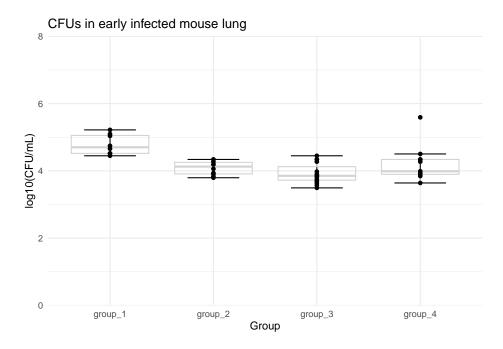

5.9 Sample ANOVA

```
cfu_stats <- tidy_cfu_meta_joined %>%
  group_by(organ) %>%
 nest() %>%
  mutate(aov_result = map(data, ~aov(CFUs_per_ml ~ groups, data = .x)),
         tukey_result = map(aov_result, TukeyHSD),
         tidy_tukey = map(tukey_result, broom::tidy)) %>%
  unnest(tidy_tukey, .drop = TRUE) %>%
  separate(contrast, into = c("contrast1", "contrast2"), sep = "-") %>%
  select(-data, -aov_result, -tukey_result, -term, -null.value)# %>%
## Warning: The `.drop` argument of `unnest()` is deprecated as of tidyr 1.0.0.
## All list-columns are now preserved.
## This warning is displayed once every 8 hours.
## Call `lifecycle::last_lifecycle_warnings()` to see where this warning was generated.
  # filter(adj.p.value <= 0.05)</pre>
cfu\_stats
## # A tibble: 9 x 7
## # Groups:
              organ [2]
     organ contrast1 contrast2 estimate conf.low conf.high adj.p.value
    <chr> <chr> <chr>
                               <dbl> <dbl>
                                                     <dbl>
```

```
## 1 lung group_2 group_1
                                              -60953. -138742.
                                                                            16836.
                                                                                             0.171
## 2 lung group_3 group_1
                                              -63903. -135699.
                                                                            7893.
                                                                                             0.0963
## 3 lung group_4 group_1 -26214. -102416.
                                                                            49987.
                                                                                             0.793
## 4 lung group_3 group_2 -2950. -69900.
                                                                          64000.
                                                                                             0.999
## 5 lung group_4 group_2 34739. -36915. 106393.

## 6 lung group_4 group_3 37689. -27410. 102787.

## 7 spleen group_2 group_1 -6565 -13529. 399.

## 8 spleen group_3 group_1 -7310 -13341. -1279.

## 9 spleen group_3 group_2 -745. -6776. 5286.
                                                                                             0.569
                                                                                             0.417
                                                                                             0.0656
                                                                                             0.0178
                                                                                             0.943
```

5.10 Save processed data to database

5.11 Example two

Enzyme-linked immunosorbest assay (ELISA)

ELISA is a standard molecular biology assay for detecting and quantifying a variety of compounds, including peptides, proteins, and antibodies in a sample. The sample could be serum, plasma, or bronchoalveolar lavage fluid (BALF).

6.0.0.1 Importance of ELISA

An antigen-specific reaction in the host results in the production of antibodies, which are proteins found in the blood. In the event of an infectious disease, it aids in the detection of antibodies in the body. ELISA is distinguishable from other antibody-assays in that it produces quantifiable findings and separates non-specific from specific interactions by serial binding to solid surfaces, which is often a polystyrene multiwell plate.

In IMPAc-TB project, it is crucial to evaluate the if the vaccine is eliciting humoral immunity and generating antibodies against vaccine antigen. ELISA will be used to determine the presence of Immunoglobulin (Ig) IgG, IgA, and IgM in the serum different time points post-vaccination.

6.0.0.2 Principle of ELISA

ELISA is based on the principle of antigen-antibody interaction. An antigen must be immobilized on a solid surface and then complexed with an enzymelinked antibody in an ELISA. The conjugated enzyme's activity is evaluated by incubating it with a substrate to yield a quantifiable result, which enables detection. There are four basic steps of ELISA:

- 1. Coating multiwell plate with antigen/antibody: This step depends on what we want to detect the sample. If we need to evaluate the the presence of antibody, the plate will be coated with the antigen, and vice versa. To coat the plate, a fixed concentration of antigen (protein) is added to a 96 well high-binding plate (charged plate). Plate is incubated over night with the antigen at 4 degree celsius (as proteins are temperature sensitive) so that antigens are completely bound to the well.
- 2. Blocking: It is possible that not each and every site of the well is coated with the targeted antigen, and there could be uncovered areas. It is important to block those empty spaces so that primary antibody (which we will add to the next step) binds to these spaces and give us false positive results. For this, microplate well surface-binding sites are blocked with an unrelated protein or other substance. Most common blocking agents are bovine serum albumin, skim milk, and casein. One of the best blocking agents is to use the serum from the organism in which your secondary (detection antibody) is raised. For example, if the secondary antibody is raised in goat, then we can use goat serum as a blocking agent.
- **3. Probing**: Probing is the step where we add sample containing antibodies that we want to detect. This will be the primary antibody. If the antibodies against the antigen (which we have coated) are present in the sample, it will bind to the antigen with high affinity.
- 4. Washing: After the incubation of sample containing primary antibody, the wells are washed so that any unbound antibody is washed away. Washing solution contains phosphate buffer saline + 0.05% tween-20 (a mild detergent). 0.05% tween-20 washes away all the non-specific interactions as those are not strong, but keeps all the specific interaction as those are strong and cannot be detached with mild detergent.
- **5. Detection**: To detect the presence of antibody-antigen complex, a secondary antibody labelled with an enzyme (usually horseradish peroxidase) is added to the wells, incubated and washed.
- **6. Signal Measurement**: Finally to detect "if" and "how much" of the antibody is present, a chromogenic substrate (like 3,3',5,5'-Tetramethylbenzidine) is added to the wells, which can be cleaved the enzyme that is tagged to the secondary antibody. The color compund is formed after the addition of the substrate, which is directly proportional to the amount of antibody present in the sample. The plate is read on a plate reader, where color is converted to numbers.

6.1 Read in data from excel file

library(readxl)
library(tidyverse)

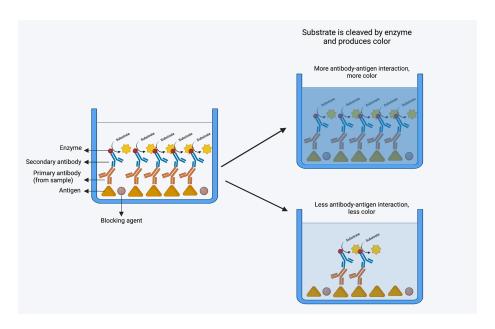


Figure 6.1: A caption

```
library(minpack.lm)
library(broom)
library(purrr)
elisa_raw_data <- read_excel("DATA/elisa_s1_07-25-20.xlsx", sheet = "S1", col_names = FALSE, rar
## New names:
## * `` -> ...1
## * `` -> ...2
## * `` -> ...3
## * `` -> ...4
## * `` -> ...5
## * ...
head(elisa_raw_data)
## # A tibble: 6 x 12
##
               ...2 ...3 ...4
    ...1
                              ...5 ...6 ...7
                                            ...8 ...9 ...10 ...11 ...12
              ## 2 7.900000000~ 0.098 0.069 6.80~ 0.115 0.202 5.89~ 0.134 0.069 0.106 0.05   0.075
## 3 8.899999999~ 0.133 0.119 OVRF~ 3.87 2.32 OVRF~ 3.85 2.12 OVRF~ 3.21 1.02
              3.46 1.16 OVRF~ 3.80 2.36
## 4 OVRFLW
                                      OVRF~ 3.70 1.49 OVRF~ 3.68 1.63
## 5 3.815999999~ 1.82 0.446 3.89~ 3.42 1.13 OVRF~ 2.33 0.608 OVRF~ 3.41 1.10
```

```
## 6 OVRFLW 3.69 1.43 OVRF~ 3.66 1.27 3.839 1.74 0.444 2.49~ 0.637 0.704
```

6.1.1 Tidy the data

```
# Convert all columns to numeric
elisa_raw_data_numeric <- elisa_raw_data %>%
  mutate_if(is.character, as.numeric)
## Warning in mask$eval_all_mutate(quo): NAs introduced by coercion
# pivot longer the data
elisa_raw_data_tidy <- pivot_longer(data = elisa_raw_data_numeric, cols = "...1":"...1
# remove "..." from the first column
elisa_raw_data_tidy$well_id <- str_replace(elisa_raw_data_tidy$well_id, "...", "")
# Add new column to the data_frame
elisa_raw_data_tidy_new <- elisa_raw_data_tidy %>%
  mutate(name = rep(LETTERS[1:8], each = 12))
elisa_raw_data_tidy_new <- elisa_raw_data_tidy_new %>%
  mutate(well_id = paste0(name, well_id)) %>%
  select(-name)
head(elisa_raw_data_tidy_new)
## # A tibble: 6 x 2
     well_id od_450nm
##
     <chr>
##
                <dbl>
## 1 A1
                0.052
## 2 A2
                0.05
## 3 A3
                0.069
## 4 A4
                0.063
## 5 A5
                0.061
## 6 A6
                0.122
```

6.1.2 Read in second data set

```
elisa_label_data <- read_excel("DATA/elisa_s1_07-25-20.xlsx", sheet = "S1", col_names = FALSE, 1
## New names:
## * `` -> ...1
## * `` -> ...2
## * `` -> ...3
## * `` -> ...4
## * `` -> ...5
## * ...
head(elisa_label_data)
## # A tibble: 6 x 12
                                           ...3 ...4 ...5 ...6 ...7 ...8 ...9 ...10 ...11 ...12
      . . . 1
                              ...2
                              <chr> <chr< <chr> <chr< <chr> <chr< <chr> <chr
        <chr>
                              secon~ naïv~ 1A-1~ 1A-1~ 1A-1~ 1A-2~ 1A-2~ 1A-2~ 1A-3~ 1A-3~ 1A-3~
## 1 blank
## 2 1A-4 (1/250 1A-4 ~ 1A-4~ 1B-1~ 1B-1~ 1B-1~ 1B-2~ 1B-2~ 1B-2~ 1B-2~ 1B-3~ 1B-3~ 1B-3~
## 3 1B-4 (1/250 1B-4 ~ 1B-4~ 2A-1~ 2A-1~ 2A-1~ 2A-2~ 2A-2~ 2A-2~ 2A-2~ 2A-3~ 2A-3~ 2A-3~
## 4 2B-1 (1/250 2B-1 ~ 2B-1~ 2B-2~ 2B-2~ 2B-2~ 2B-3~ 2B-3~ 2B-3~ 2B-4~ 2B-4~ 2B-4~
## 5 3A-1 (1/250 3A-1 ~ 3A-1~ 3A-2~ 3A-2~ 3A-2~ 3A-3~ 3A-3~ 3A-3~ 3A-4~ 3A-4~ 3A-4~
## 6 3B-1 (1/250 3B-1 ~ 3B-1~ 3B-2~ 3B-2~ 3B-2~ 3B-3~ 3B-3~ 3B-3~ 3B-4~ 3B-4~ 3B-4~
# pivot longer the data
elisa_label_data_tidy <- pivot_longer(data = elisa_label_data, cols = "...1":"...12", names_to =
# remove "..." from the first column
elisa_label_data_tidy$well_id <- str_replace(elisa_label_data_tidy$well_id, "...", "")
# Add new column to the data_frame
elisa_label_data_tidy_new <- elisa_label_data_tidy %>%
   mutate(name = rep(LETTERS[1:8], each = 12))
elisa_label_data_tidy_new <- elisa_label_data_tidy_new %>%
   mutate(well_id = paste0(name, well_id)) %>%
   select(-name)
head(elisa_label_data_tidy_new)
## # A tibble: 6 x 2
## well id information
        <chr> <chr>
## 1 A1
                       blank
```

```
## 2 A2 secondary
## 3 A3 naïve (1/250)
## 4 A4 1A-1 (1/250
## 5 A5 1A-1 (1/1250
## 6 A6 1A-1 (1/6250
```

6.2 #Join the OD table with the information table

```
elisa_data = elisa_raw_data_tidy_new %>% inner_join(elisa_label_data_tidy_new, by="well
head(elisa_data)
## # A tibble: 6 x 3
    well_id od_450nm information
               <dbl> <chr>
##
## 1 A1
               0.052 blank
## 2 A2
               0.05 secondary
## 3 A3
               0.069 naïve (1/250)
## 4 A4
               0.063 1A-1 (1/250
## 5 A5
               0.061 1A-1 (1/1250
## 6 A6
               0.122 1A-1 (1/6250
```

6.2.1 Separate the information table into sample ID and dilution columns

```
tidy_elisa_data <- separate(elisa_data, col = "information", into = c("sample_id", "di
## Warning: Expected 2 pieces. Missing pieces filled with `NA` in 2 rows [1, 2].
head(tidy_elisa_data)
## # A tibble: 6 x 4
   well_id od_450nm sample_id dilution
    <chr>
               <dbl> <chr>
##
                                 <chr>
## 1 A1
               0.052 "blank"
                                 <NA>
## 2 A2
               0.05 "secondary" <NA>
               0.069 "naïve "
## 3 A3
                                 1/250)
## 4 A4
               0.063 "1A-1 "
                                 1/250
## 5 A5
               0.061 "1A-1 "
                                 1/1250
               0.122 "1A-1 "
## 6 A6
                                 1/6250
tidy_elisa_data <- tidy_elisa_data %>%
 mutate(dilution = str_extract(dilution, "(/)[0-9]+"),
        dilution = str_replace(dilution, "/", ""),
```

```
dilution = as.numeric(dilution))

tidy_elisa_data <- tidy_elisa_data %>%
   select(well_id, sample_id, dilution, od_450nm)

head(tidy_elisa_data)
```

```
## # A tibble: 6 x 4
##
     well_id sample_id
                           dilution od_450nm
##
     <chr>
              <chr>
                              <dbl>
                                        <dbl>
## 1 A1
              "blank"
                                 NA
                                        0.052
                                        0.05
## 2 A2
              "secondary"
                                 NA
              "naïve "
## 3 A3
                                250
                                        0.069
## 4 A4
              "1A-1 "
                                250
                                        0.063
              "1A-1 "
## 5 A5
                               1250
                                        0.061
              "1A-1 "
## 6 A6
                               6250
                                        0.122
```

6.3 ELISA data analysis optimization

ELISA data can be analyzed in different ways based on how the data is acquired. There are a a few examples of the type of ELISA data:

- 1. With standard curve: ELISA can be used to determine the concentrations of the antigen and antibody. This type of ELISA data usually have a standard curve with different concentrations of the known analyte and the concentration in the sample is determined by extrapolating the unknown values in the curve. This type of assay is straightforward, easy to interpret and are more robust.
- 2. Without standard curve: Usually vaccine studies involve investigating the presence of high-affinity (and novel) antibodies against the vaccine antigens. Therefore, plotting a standard curve is not feasible as there is no previous information available for antibody concentration or type of antibody. Also, because antibody response to a vaccine will differ depending on the individual, it is not practical to generate a calibration curve from which absolute concentrations can be extrapolated. For this type of ELISA, quantification of the antibody titers is performed using serial dilutions of the test samples, and analysis can be performed using the following three methods:
 - 1. Fitting sigmoid model
 - 2. Endpoint titer method 3: Absorbance summation method


Let's have a look at these methods, how we can apply these methods in our data, and R-based packages that we can utilize to perform tis analysis.

6.3.1 1. Fitting sigmoid model:

In this model, we will perform a 8-10 point serial dilution of our sample and plot a 4 parameter sigmoidal curve. Example data:

```
elisa_example_data <- read_excel("DATA/example_elisa_data.xlsx")</pre>
# separate 1/
elisa_example_data <- separate(elisa_example_data, col = "dilution", into = c("numerat
elisa_example_data <- elisa_example_data %>%
  mutate_if(is.character, as.numeric)
elisa_example_data$dilution <- elisa_example_data$numerator / elisa_example_data$denom
elisa_example_data <- elisa_example_data %>%
  mutate(log_dilution = log(dilution, base = 3))
head(elisa_example_data)
## # A tibble: 6 x 5
     numerator denomenator absorbance dilution log_dilution
##
         <dbl>
                     <dbl>
                                <dbl>
                                          <dbl>
                                                       <dbl>
## 1
                                      0.0333
                                                       -3.10
             1
                        30
                                  4
## 2
             1
                        90
                                  3.73 0.0111
                                                       -4.10
## 3
             1
                       270
                                  2.34 0.00370
                                                       -5.10
## 4
                       810
                                 1.1 0.00123
                                                       -6.10
             1
## 5
                      2430
                                  0.51 0.000412
                                                       -7.10
             1
## 6
                      7290
                                  0.22 0.000137
                                                       -8.10
mod_1 <- nlsLM(absorbance ~ ((a-d)/(1+(log_dilution/c)^b)) + d,</pre>
data = elisa_example_data,
start = list (a = 4, d = 0, c = -5, b = 1))
\# a = maximum \ absorbance
\# d = minimum \ absobance
# c = point of maximum growith (dilition where the curve is straight line)
# b = slope at c
mod_1
## Nonlinear regression model
     model: absorbance \sim ((a - d)/(1 + (\log_dilution/c)^b)) + d
##
      data: elisa_example_data
##
## 4.12406 0.04532 -5.31056 7.62972
```

```
## residual sum-of-squares: 0.02221
##
## Number of iterations to convergence: 9
## Achieved convergence tolerance: 1.49e-08
summary(mod_1)
##
## Formula: absorbance \sim ((a - d)/(1 + (\log_dilution/c)^b)) + d
## Parameters:
## Estimate Std. Error t value Pr(>|t|)
## a 4.12406 0.05820 70.860 1.75e-12 ***
## d 0.04532 0.02268 1.998 0.0808.
## c -5.31056 0.03933 -135.037 1.01e-14 ***
## b 7.62972 0.35854 21.280 2.50e-08 ***
## ---
## Signif. codes: 0 '***' 0.001 '**' 0.05 '.' 0.1 ' ' 1
##
## Residual standard error: 0.05269 on 8 degrees of freedom
##
## Number of iterations to convergence: 9
## Achieved convergence tolerance: 1.49e-08
6311 Fitted model curve
tidy_params <- mod_1 %>% tidy()
a <- tidy_params$estimate[tidy_params$term == "a"]</pre>
b <- tidy_params$estimate[tidy_params$term == "b"]</pre>
c <- tidy_params$estimate[tidy_params$term == "c"]</pre>
d <- tidy_params$estimate[tidy_params$term == "d"]</pre>
elisa_example_data <- elisa_example_data %>%
 mutate(fitted = predict(mod_1))
elisa_example_data <- elisa_example_data %>%
  mutate(fitted = predict(mod_1))
elisa_example_data %>%
  ggplot(aes(x = log_dilution, y = absorbance)) +
  geom point() +
 geom_line(aes(y=fitted), color = "blue")
```

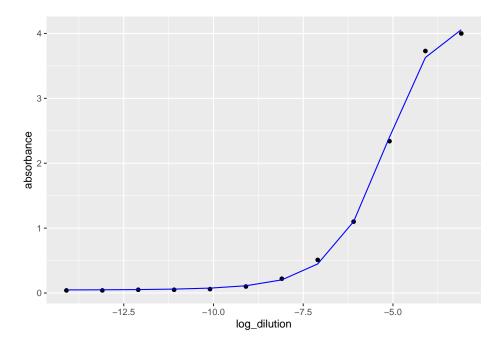

6.3.2 2. Endpoint titer method

The endpoint titer approach chooses an absorbance value just above the background noise (or the lower asymptotic level). The highest dilution with an absorbance greater than this predetermined value is the endpoint titer. This method is based on the assumption that a sample with a higher protein concentration will require a higher dilution factor to achieve an absorbance just above the level of background noise.

```
endpoint_titer <- c * (((a - d) / (0.2 - d)) - 1) ^ (1 / b)
summary(endpoint_titer)

## Min. 1st Qu. Median Mean 3rd Qu. Max.
## -8.113 -8.113 -8.113 -8.113 -8.113
endpoint_titer</pre>
```

[1] -8.113285

6.3.3 Absorbance summation method

In this model of data analysis, we sum all the absorbance values from each sample to obtain one value. This value is termed as absorption summation (AS). Using the above data, the AS will be calculated as below:

```
AS = 0.04 + 0.04 + 0.05 + 0.05 + 0.06 + 0.1 + 0.22 + 0.51 + 1.1 + 2.34 + 3.73 + 4.0

AS = ## [1] 12.24
```

6.4 Use the the data analysis models in our data

The presented data is from a mouse study. In this data, presence of IgG antibody has been evaluated against receptor binding domain (RBD) of SARS-CoV-2 virus in two different groups of mice. We need to elucidate which group has higher concentration of the antibodies.

6.4.1 read in the data

```
elisa_data <- read_excel("DATA/elisa_data_serial_dilution.xlsx")</pre>
```

denomenator = as.numeric(denomenator))

6.4.2 make the tidy data

```
elisa_data <- pivot_longer(data = elisa_data, cols = "Mouse_1": "Mouse_5", names_to = "mouse_id",
elisa_data
## # A tibble: 100 x 4
     Groups Dilution mouse_id absorbance
##
      <chr>
              <chr>>
                       <chr>
                                     <dbl>
## 1 Group 1 1/50
                       Mouse_1
                                       4.1
## 2 Group 1 1/50
                      Mouse_2
                                       3.9
## 3 Group 1 1/50
                       Mouse_3
                                       4.3
## 4 Group 1 1/50
                       Mouse_4
                                       4.2
## 5 Group 1 1/50
                       Mouse_5
                                       4
## 6 Group 1 1/100
                       Mouse_1
                                       3.9
## 7 Group 1 1/100
                       Mouse 2
                                       3.8
## 8 Group 1 1/100
                       Mouse_3
                                       3.6
## 9 Group 1 1/100
                       Mouse_4
                                       3.7
## 10 Group 1 1/100
                                       3.8
                       Mouse_5
## # ... with 90 more rows
# separate dilution column and convert it to log2
elisa data <- separate(elisa data, col = "Dilution", into = c("numerator", "denomenator"), sep =
elisa_data <- elisa_data %>%
  transform(numerator = as.numeric(numerator),
```

```
elisa_data <- elisa_data %>%
  mutate(dilution = elisa_data$numerator / elisa_data$denomenator)
elisa_data <- elisa_data %>%
  mutate(log_dilution = log2(dilution))
elisa_data
```

```
##
        Groups numerator denomenator mouse_id absorbance
                                                                dilution log_dilution
## 1
       Group 1
                                    50
                                        Mouse 1
                                                       4.10 2.00000e-02
                                                                             -5.643856
                        1
## 2
       Group 1
                                    50
                                        Mouse_2
                                                       3.90 2.00000e-02
                                                                             -5.643856
                        1
## 3
       Group 1
                        1
                                    50
                                        Mouse_3
                                                       4.30 2.00000e-02
                                                                             -5.643856
## 4
                                    50
       Group 1
                                        Mouse_4
                                                       4.20 2.00000e-02
                                                                             -5.643856
                        1
## 5
       Group 1
                                    50
                                        Mouse_5
                                                       4.00 2.00000e-02
                                                                             -5.643856
                        1
## 6
       Group 1
                                   100
                                        Mouse_1
                                                       3.90 1.00000e-02
                                                                             -6.643856
                        1
## 7
       Group 1
                        1
                                   100
                                        Mouse 2
                                                       3.80 1.00000e-02
                                                                             -6.643856
                                   100
## 8
       Group 1
                        1
                                        Mouse_3
                                                       3.60 1.00000e-02
                                                                             -6.643856
## 9
       Group 1
                        1
                                   100
                                        Mouse_4
                                                       3.70 1.00000e-02
                                                                             -6.643856
## 10
                                   100
                                                       3.80 1.00000e-02
       Group 1
                        1
                                        Mouse_5
                                                                             -6.643856
## 11
       Group 1
                        1
                                   200
                                        Mouse 1
                                                       3.42 5.00000e-03
                                                                             -7.643856
## 12
       Group 1
                        1
                                   200
                                        Mouse 2
                                                       3.35 5.00000e-03
                                                                            -7.643856
## 13
       Group 1
                        1
                                   200
                                        Mouse 3
                                                       3.56 5.00000e-03
                                                                             -7.643856
                                   200
## 14
       Group 1
                        1
                                        Mouse 4
                                                       3.28 5.00000e-03
                                                                             -7.643856
## 15
       Group 1
                        1
                                   200
                                        Mouse 5
                                                       3.29 5.00000e-03
                                                                             -7.643856
## 16
       Group 1
                        1
                                   400
                                        Mouse_1
                                                       2.89 2.50000e-03
                                                                             -8.643856
                                                                             -8.643856
## 17
       Group 1
                        1
                                   400
                                        Mouse 2
                                                       2.98 2.50000e-03
## 18
       Group 1
                        1
                                   400
                                        Mouse 3
                                                       2.64 2.50000e-03
                                                                             -8.643856
                                   400
                                        Mouse 4
## 19
       Group 1
                        1
                                                       2.78 2.50000e-03
                                                                            -8.643856
## 20
       Group 1
                        1
                                   400
                                        Mouse_5
                                                       2.91 2.50000e-03
                                                                            -8.643856
                                                                             -9.643856
## 21
       Group 1
                                   800
                                        Mouse_1
                                                       2.12 1.25000e-03
                        1
## 22
       Group 1
                        1
                                   800
                                        Mouse_2
                                                       2.18 1.25000e-03
                                                                             -9.643856
## 23
                        1
                                   800
       Group 1
                                        Mouse_3
                                                       2.17 1.25000e-03
                                                                             -9.643856
                        1
                                   800
                                                       2.14 1.25000e-03
## 24
       Group 1
                                        Mouse 4
                                                                             -9.643856
## 25
       Group 1
                                   800
                                        Mouse_5
                                                       2.21 1.25000e-03
                                                                            -9.643856
                        1
## 26
       Group 1
                        1
                                  1600
                                        Mouse_1
                                                       1.89 6.25000e-04
                                                                           -10.643856
## 27
       Group 1
                        1
                                  1600
                                        Mouse_2
                                                       1.98 6.25000e-04
                                                                           -10.643856
## 28
       Group 1
                        1
                                  1600
                                        Mouse 3
                                                       1.76 6.25000e-04
                                                                           -10.643856
## 29
       Group 1
                                        Mouse 4
                                                       1.69 6.25000e-04
                                                                           -10.643856
                        1
                                  1600
## 30
       Group 1
                        1
                                  1600
                                        Mouse 5
                                                       1.82 6.25000e-04
                                                                           -10.643856
       Group 1
                        1
## 31
                                  3200
                                        Mouse_1
                                                       0.97 3.12500e-04
                                                                           -11.643856
                                        Mouse_2
## 32
       Group 1
                        1
                                  3200
                                                       0.85 3.12500e-04
                                                                           -11.643856
## 33
       Group 1
                        1
                                  3200
                                        Mouse_3
                                                       0.88 3.12500e-04
                                                                           -11.643856
## 34
       Group 1
                                  3200
                                        Mouse 4
                                                       0.84 3.12500e-04
                        1
                                                                           -11.643856
## 35
       Group 1
                                  3200
                                        Mouse 5
                                                       0.91 3.12500e-04
                                                                           -11.643856
                                                                           -12.643856
## 36
       Group 1
                        1
                                  6400
                                        Mouse 1
                                                       0.52 1.56250e-04
```

##	37	Group	1	1	6400	Mouse 2	0.48 1.56250e	e-04 -12.643856
##		Group		1	6400	Mouse_3	0.42 1.562506	
	39	Group		1	6400	Mouse_4	0.48 1.56250e	
	40	Group		1	6400	Mouse 5	0.47 1.56250e	
	41	Group		1	12800	Mouse_1	0.26 7.81250	
	42	Group		1	12800	Mouse_1	0.20 7.81250e	
	43	Group		1	12800	Mouse_3	0.27 7.81250e	
	44	-		1	12800	_	0.26 7.81250e	
	45	Group		1	12800	Mouse_4	0.28 7.81250e	
	46	Group		1	25600	Mouse_5 Mouse_1	0.12 3.90625e	
	47	Group Group		1	25600		0.12 3.90025e	
	48	Group		1	25600	Mouse_2	0.10 3.90025e	
##		-		1	25600	Mouse_3 Mouse_4	0.16 3.906256	
##		Group		1	25600	Mouse_5	0.14 3.906256	
##		Group		1	50	_	3.20 2.00000e	
##		Group		1	50	Mouse_1	3.40 2.00000e	
##		Group Group		1	50	Mouse_2	3.60 2.00000e	
##		Group		1	50	Mouse_3	3.20 2.00000e	
##		Group		1	50	Mouse_4	3.70 2.00000e	
##		-		1	100	Mouse_5	2.90 1.00000e	
##		Group Group		1	100	Mouse_1 Mouse 2	2.80 1.00000e	
	58	Group		1	100	Mouse_3	2.70 1.00000e	
	59	Group		1	100	Mouse_4	2.80 1.00000e	
##		Group		1	100	Mouse_5	2.80 1.00000e	
	61	Group		1	200	Mouse_1	2.32 5.000006	
	62	Group		1	200	Mouse_2	2.35 5.00000e	
##		Group		1	200	Mouse_3	2.56 5.00000e	
##		Group		1	200	Mouse_4	2.28 5.000006	
##		Group		1	200	Mouse_5	2.29 5.000006	
	66	Group		1	400	Mouse_1	1.89 2.50000e	
	67	Group		1	400	Mouse_2	1.98 2.50000e	
	68	Group		1	400	Mouse_3	1.64 2.50000e	
	69	Group		1	400	Mouse_4	1.78 2.50000e	
	70	Group		1	400	Mouse_5	1.91 2.50000e	
	71	Group		1	800	Mouse_1	1.12 1.25000e	
	72	Group		1	800	Mouse_2	1.18 1.25000e	
##	73	Group		1	800	Mouse_3	1.17 1.25000e	
	74	Group	_	1	800	Mouse_4	1.14 1.25000e	
##	75	Group	2	1	800	Mouse_5	1.21 1.25000e	-03 -9.643856
##	76	Group		1	1600	Mouse_1	0.89 6.25000e	-04 -10.643856
##	77	Group		1	1600	Mouse_2	0.98 6.25000e	-04 -10.643856
##	78	Group	2	1	1600	Mouse_3	0.76 6.25000€	e-04 -10.643856
##	79	Group		1	1600	Mouse_4	0.69 6.25000e	e-04 -10.643856
##	80	Group	2	1	1600	Mouse_5	0.82 6.25000e	e-04 -10.643856
##	81	Group	2	1	3200	Mouse_1	0.47 3.12500€	e-04 -11.643856
##	82	${\tt Group}$	2	1	3200	Mouse_2	0.45 3.12500e	e-04 -11.643856

```
## 83
      Group 2
                               3200
                                    Mouse_3
                                                  0.48 3.12500e-04
                                                                     -11.643856
                      1
                                                  0.44 3.12500e-04
## 84
      Group 2
                      1
                               3200
                                    Mouse_4
                                                                     -11.643856
## 85
      Group 2
                               3200
                                    Mouse_5
                      1
                                                  0.41 3.12500e-04
                                                                     -11.643856
## 86
      Group 2
                      1
                               6400
                                    Mouse 1
                                                  0.20 1.56250e-04
                                                                     -12.643856
## 87
      Group 2
                      1
                               6400
                                    Mouse 2
                                                  0.28 1.56250e-04
                                                                     -12.643856
## 88
      Group 2
                      1
                               6400
                                    Mouse 3
                                                  0.22 1.56250e-04
                                                                    -12.643856
## 89
      Group 2
                               6400
                                    Mouse 4
                                                  0.28 1.56250e-04
                     1
                                                                    -12.643856
## 90
      Group 2
                               6400
                                    Mouse_5
                                                 0.27 1.56250e-04
                                                                    -12.643856
                      1
                                                 0.16 7.81250e-05
## 91
      Group 2
                      1
                              12800
                                    Mouse_1
                                                                     -13.643856
                                                                    -13.643856
## 92
      Group 2
                      1
                              12800
                                    Mouse 2
                                                  0.11 7.81250e-05
## 93
      Group 2
                     1
                              12800
                                    Mouse 3
                                                  0.17 7.81250e-05 -13.643856
## 94
      Group 2
                      1
                              12800
                                    Mouse 4
                                                  0.16 7.81250e-05
                                                                    -13.643856
## 95
      Group 2
                      1
                              12800
                                    Mouse_5
                                                  0.18 7.81250e-05
                                                                     -13.643856
                                                  0.02 3.90625e-05
                                    Mouse_1
## 96
      Group 2
                     1
                              25600
                                                                    -14.643856
## 97
      Group 2
                     1
                              25600
                                    Mouse_2
                                                  0.06 3.90625e-05
                                                                     -14.643856
                                                                     -14.643856
## 98
      Group 2
                      1
                              25600
                                    Mouse_3
                                                  0.04 3.90625e-05
## 99
      Group 2
                              25600
                                    Mouse 4
                                                  0.08 3.90625e-05
                                                                     -14.643856
                      1
## 100 Group 2
                      1
                              25600
                                    Mouse_5
                                                  0.06 3.90625e-05
                                                                     -14.643856
# converting data into dataframe
elisa_data_df <- elisa_data %>%
  group_by(Groups, mouse_id) %>%
  summarize(log_dilution = log_dilution, absorbance = absorbance)
## `summarise()` has grouped output by 'Groups', 'mouse_id'. You can override using
## the `.groups` argument.
elisa_data_df
## # A tibble: 100 x 4
              Groups, mouse_id [10]
## # Groups:
     Groups mouse_id log_dilution absorbance
##
      <chr>
             <chr>
                                        <dbl>
                            <dbl>
## 1 Group 1 Mouse_1
                            -5.64
                                         4.1
                            -6.64
## 2 Group 1 Mouse 1
                                        3.9
                            -7.64
## 3 Group 1 Mouse 1
                                        3.42
## 4 Group 1 Mouse_1
                            -8.64
                                        2.89
## 5 Group 1 Mouse_1
                            -9.64
                                        2.12
## 6 Group 1 Mouse_1
                            -10.6
                                        1.89
## 7 Group 1 Mouse_1
                            -11.6
                                        0.97
## 8 Group 1 Mouse 1
                            -12.6
                                        0.52
```

0.26

0.12

-13.6

-14.6

```
elisa_data_nested <- elisa_data %>%
  group_by(Groups, mouse_id) %>%
```


9 Group 1 Mouse_1

10 Group 1 Mouse_1

... with 90 more rows

```
nest()
elisa_data_nested
```

```
## # A tibble: 10 x 3
## # Groups: Groups, mouse_id [10]
##
     Groups mouse_id data
##
     <chr>
             <chr> <chr>
## 1 Group 1 Mouse_1 <tibble [10 x 5]>
## 2 Group 1 Mouse_2 <tibble [10 x 5]>
## 3 Group 1 Mouse_3 <tibble [10 x 5]>
## 4 Group 1 Mouse_4 <tibble [10 x 5]>
## 5 Group 1 Mouse_5 <tibble [10 x 5]>
## 6 Group 2 Mouse_1 <tibble [10 x 5]>
## 7 Group 2 Mouse_2 <tibble [10 x 5]>
## 8 Group 2 Mouse_3 <tibble [10 x 5]>
## 9 Group 2 Mouse_4 <tibble [10 x 5]>
## 10 Group 2 Mouse_5 <tibble [10 x 5]>
# plot the curves to evaluate the a, d, c, and b
elisa_data %>%
 ggplot(aes(x = log_dilution, y = absorbance)) +
 geom_point() +
 geom_line() +
 facet_wrap(Groups ~ mouse_id)
```

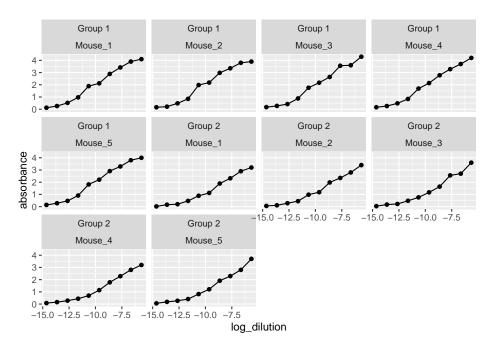

Based on the curve, the values are:

$$a = 4, d = 0 c = 2 b = 1$$

6.5 Creating a function for fitting model

```
# Creating a function for fitted model

fitted_model_elisa <- function(df_elisa, start_a, start_d, start_c, start_b) {
  mod_1 <- nlsLM(absorbance ~ ((a-d)/(1+(log_dilution/c)^b)) + d,
  data = df_elisa,
  start = list(a = start_a, d = start_d, c = start_c, b = start_b))
  return(mod_1)
}</pre>
```

6.5.1 Fitting the model into the data

```
fitted_model_elisa(elisa_data_nested$data[[1]], start_a = 4, start_d = 0, start_c = -8
## Nonlinear regression model
## model: absorbance ~ ((a - d)/(1 + (log_dilution/c)^b)) + d
## data: df_elisa
## a d c b
## 4.3070 -0.6009 -10.2577 5.2893
```

```
## residual sum-of-squares: 0.1199
##
## Number of iterations to convergence: 7
## Achieved convergence tolerance: 1.49e-08
```

6.5.2 Apply the fitted model function to the whole dataframe

```
elisa_fitted_data <- elisa_data_nested %>%
 mutate(fitted_data = purrr::map(data, ~ fitted_model_elisa(.x, start_a = 4, start_d = 0, start_
elisa_fitted_data
## # A tibble: 10 x 4
## # Groups:
              Groups, mouse_id [10]
##
     Groups mouse_id data
                                        fitted_data
##
     <chr>
             <chr>
                      st>
                                        t>
## 1 Group 1 Mouse_1 <tibble [10 x 5]> <nls>
## 2 Group 1 Mouse_2 <tibble [10 x 5]> <nls>
## 3 Group 1 Mouse_3 <tibble [10 x 5]> <nls>
## 4 Group 1 Mouse_4 <tibble [10 x 5]> <nls>
## 5 Group 1 Mouse_5 <tibble [10 x 5]> <nls>
## 6 Group 2 Mouse_1 <tibble [10 x 5]> <nls>
## 7 Group 2 Mouse_2 <tibble [10 x 5]> <nls>
## 8 Group 2 Mouse_3 <tibble [10 x 5]> <nls>
## 9 Group 2 Mouse_4 <tibble [10 x 5]> <nls>
## 10 Group 2 Mouse_5 <tibble [10 x 5]> <nls>
```

6.5.3 Take out the summary of the data

```
elisa_fitted_data_summary <- elisa_fitted_data %>%
   mutate(fitted_data_summary = purrr::map(fitted_data, broom::glance))
#unnested_elisa_fitted_data <- elisa_fitted_data_summary %>%
# unnest(elisa_fitted_data_summary) %>%
# ungroup() %>%
# dplyr::select(Groups, mouse_id, fitted_data)
# unnested_elisa_fitted_data$fitted_data[[1]]
# Another way
elisa_fitted_data_summary %>%
   unnest(fitted_data_summary)
```

```
## # A tibble: 10 x 13
## # Groups:
              Groups, mouse_id [10]
##
     Groups mouse_id data
                              fitted_data sigma isConv
                                                           finTol logLik
                                                                            AIC
                                          <dbl> <lgl>
             <chr>
                              t>
                                                            <dbl> <dbl>
##
     <chr>
                     t>
                                                                          <dbl>
                                         0.141 TRUE
## 1 Group 1 Mouse_1 <tibble> <nls>
                                                          1.49e-8
                                                                   7.93 - 5.86
##
  2 Group 1 Mouse 2 <tibble> <nls>
                                         0.188 TRUE
                                                          1.49e-8 5.10 -0.208
   3 Group 1 Mouse_3 <tibble> <nls>
                                         0.204 TRUE
                                                          1.49e-8 4.25
                                                                         1.49
   4 Group 1 Mouse_4 <tibble> <nls>
##
                                         0.132 TRUE
                                                          1.49e-8 8.65 -7.30
   5 Group 1 Mouse 5 <tibble> <nls>
                                         0.134 TRUE
                                                          1.49e-8
                                                                 8.50 -6.99
##
   6 Group 2 Mouse_1 <tibble> <nls>
##
                                         0.0838 TRUE
                                                          1.49e-8 13.2 -16.3
                                                                  9.71 -9.41
  7 Group 2 Mouse 2 <tibble> <nls>
                                         0.118 TRUE
                                                          1.49e-8
## 8 Group 2 Mouse_3 <tibble> <nls>
                                         0.154 TRUE
                                                          1.49e-8
                                                                  7.07 - 4.14
  9 Group 2 Mouse_4 <tibble> <nls>
                                         0.0430 TRUE
                                                          1.49e-8 19.8 -29.7
## 10 Group 2 Mouse_5 <tibble> <nls>
                                         0.126 TRUE
                                                          1.49e-8 9.11 -8.23
## # ... with 4 more variables: BIC <dbl>, deviance <dbl>, df.residual <int>,
      nobs <int>
```

6.5.4 Creating a function for endpoint titer model and fitting the model into the data

Flow cytometry