# ImpactTB/BAA: Standard Operating Procedures for Data Analysis

Colorado State University Coding Team

2022-07-07

# Contents

1	Overview	5
<b>2</b>	Introduction	7
	2.1 About the project: Immune Mechanisms of Protection against Mycobacterium tuberculosis (IMPAc-TB)	7
3	Initial mouse characteristics	9
4	Mouse Weights	11
	4.1 Read in data	12
5	Data is stored in one excel sheet, each week is one sheet named as the date -> return vector for each sheet name	13
6	can also use rio to read in the data, more streamlined	<b>15</b>
	6.1 Clean data	15
	6.2 Summary statistics	15
7	Body weight over time graph and statistics	17
8	Weight loss over time graph and statistics	19
9	Weight vs CFU	<b>21</b>
10	Weight vs ELISA results	23
11	Weight vs lesion burden	<b>25</b>

4 CONTENTS

<b>12</b>	Colo	ony forming units to determine bacterial counts	<b>27</b>				
	12.1	Data description	27				
	12.2	Read in data	28				
	12.3	Exploratory analysis and quality checks	29				
	12.4 Exploratory analysis						
	12.5	Identify a good dilution for each sample $\ldots \ldots \ldots$	30				
	12.6 Calculate CFUs from best dilution/Estimate bacterial load for each sample based on good dilution						
	12.7	Create initial report information for these data $\ldots \ldots$	31				
	12.8	Sample ANOVA	31				
	12.9	Save processed data to database	32				
	12.10Example one						
	12.11Example two						
13	Enz	yme-linked immunosorbest assay (ELISA)	33				
	13.1	Read in data from excel file	35				
	13.2	Tidy the data	36				
	13.3	Read in second data set	37				
	13.4	Join the OD table with the information table $\ \ldots \ \ldots \ \ldots$	38				
	13.5	Separate the information tabke into sample ID and dilution columns $$	39				
	13.6	ELISA data analysis optimization	40				

# Overview

Here, we have built a comprehensive guide to wet lab data collection, sample processing, and computational tool creation for robust and efficient data analysis and dissemination.

# Introduction

#### 2.1 About the project: Immune Mechanisms of Protection against Mycobacterium tuberculosis (IMPAc-TB)

The objective of the IMPAc-TB program is to get a thorough understanding of the immune responses necessary to avoid initial infection with Mycobacterium tuberculosis (Mtb), formation of latent infection, and progression to active TB illness. To achieve these goals, the National Institute of Allergy and Infectious Diseases awarded substantial funding and established multidisciplinary research teams that will analyze immune responses against Mtb in animal models (mice, guinea pigs, and non-human primates) and humans, as well as immune responses elicited by promising vaccine candidates. The contract awards establish and give up to seven years of assistance for IMPAc-TB Centers to explain the immune responses required for Mtb infection protection.

The seven centers that are part of the study are (in alphabetical order):

- 1. Colorado State University
- 2. Harvard T.H. Chan School of Public Health
- 3. Seattle Children Hospital

Colorado State University Team and role of each member: Dr. Marcela Henao-Tamayo: Principal Investigator Dr. Brendan Podell: Principal Investigator Dr. Andres Obregon-Henao: Research Scientist-III Dr. Taru S. Dutt: Research Scientist-I

# Initial mouse characteristics

Here is a review of existing methods.

# Mouse Weights

```
For example, Baazim et al. (2022) said ... (Baazim et al., 2022)
```

Mice are weighed in grams weekly to monitor clinical status as TB patients frequently display weight loss as clinical symptom associated with disease progression.

Weights are recorded in an excel worksheet.

Column titles are as follows: who\_collected date\_collected sex dob notch\_id mouse\_number weight unit cage\_number group notes

Groups included are: bcg, saline, bcg+id93, saline+id93, saline+noMtb

The notes column contains information regarding clinical observations.

good reference: https://elifesciences.org/articles/74419#s4

```
library(readxl)
## Warning: package 'readxl' was built under R version 4.1.2
library(tidyverse)
                                    ----- tidyverse 1.3.1 --
## -- Attaching packages -----
## v ggplot2 3.3.5
                    v purrr
                             0.3.4
## v tibble 3.1.6
                    v dplyr
                            1.0.7
## v tidyr 1.1.4
                    v stringr 1.4.0
## v readr
           2.1.2
                    v forcats 0.5.1
```

## Warning: package 'readr' was built under R version 4.1.2

```
## -- Conflicts ------ tidyverse_conflicts() --
## x dplyr::filter() masks stats::filter()
## x dplyr::lag() masks stats::lag()
```

#### 4.1 Read in data

Data is stored in one excel sheet, each week is one sheet named as the date -> return vector for each sheet name 14CHAPTER 5. DATA IS STORED IN ONE EXCEL SHEET, EACH WEEK IS ONE SHEET NAMED

# can also use rio to read in the data, more streamlined

#### 6.1 Clean data

#### 6.2 Summary statistics

16CHAPTER 6. CAN ALSO USE RIO TO READ IN THE DATA, MORE STREAMLINED

Body weight over time graph and statistics

#### 18CHAPTER 7. BODY WEIGHT OVER TIME GRAPH AND STATISTICS

Weight loss over time graph and statistics

#### 20 CHAPTER 8. WEIGHT LOSS OVER TIME GRAPH AND STATISTICS

Weight vs CFU

Weight vs ELISA results

Weight vs lesion burden

# Colony forming units to determine bacterial counts

#### 12.1 Data description

The data are collected in a spreadsheet with multiple sheets. The first sheet (named "[x]") is used to record some metadata for the experiment, while the following sheets are used to record CFUs counts from the plates used for samples from each organ, with one sheet per organ. For example, if you plated data from both the lung and spleen, there would be three sheets in the file: one with the metadata, one with the plate counts for the lung, and one with the plate counts for the spleen.

The metadata sheet is used to record information about the overall process of plating the data. Values from this sheet will be used in calculating the bacterial load in the original sample based on the CFU counts. This spreadsheet includes the following columns:

- organ: Include one row for each organ that was plated in the experiment. You should name the organ all in lowercase (e.g., "lung", "spleen"). You should use the same name to also name the sheet that records data for that organ for example, if you have rows in the metadata sheet for "lung" and "spleen", then you should have two other sheets in the file, one sheet named "lung" and one named "spleen", which you'll use to store the plate counts for each of those organs.
- prop\_resuspended: In this column, give the proportion of that organ that
  was plated. For example, if you plated half the lung, then in the "lung"
  row of this spread sheet, you should put 0.5 in the prop\_resuspended
  column.

- total\_resuspended\_uL: This column contains an original volume of tissue homogenate. For example, raw lung tissue is homogenized in 500 uL of PBS in a tube containing metal beads.
- og\_aliquot\_uL: 100 uL of th total\_resuspended slurry would be considered an original aliquot and is used to perform serial dilutions.
- dilution\_factor: Amount of the original stock solution that is present in the total solution, after dilution(s)
- plated\_uL: Amount of suspension + diluent plated on section of solid agar

#### 12.2 Read in data

```
library(readxl)
library(dplyr)
library(purrr)
library(tidyr)
library(stringr)
#Replace w/ path to CFU sheet
path <- c("DATA/Copy of baa_cfu_sheet.xlsx")</pre>
sheet_names <- excel_sheets(path)</pre>
sheet_names <- sheet_names[!sheet_names %in% c("metadata")]</pre>
merged_data <- list()</pre>
for(i in 1:length(sheet_names)){
  data <- read_excel(path, sheet = sheet_names[i]) %>%
    mutate(organ = paste0(sheet_names[i]))
  data <- data %>%
    #mutate(missing_col = NA) %>%
    mutate_if(is.double, as.numeric) %>%
    mutate_if(is.numeric, as.character) %>%
    pivot_longer(starts_with("dil_"), names_to = "dilution",
                  values_to = "CFUs") %>%
    mutate(dilution = str_extract(dilution, "[0-9]+"),
           dilution = as.numeric(dilution))
  merged_data[[i]] <- data
```

```
}
all_data <- bind_rows(merged_data, .id = "column_label") %>%
    select(-column_label)
```

#### 12.3 Exploratory analysis and quality checks

#### 12.4 Exploratory analysis

#### Dimensions of input data:

Based on the input data, data were collected for the following organ or organs:

The following number of mice were included for each:

The following number of replicates were recorded at each count date for each experimental group:

The following number of dilutions and dilution level were recorded for each organ:

#### People who plated and collected the data. Date or dates of counting:

Based on the input data, the plates included in these data were counted by the following person or persons: Based on the input data, the plates included in these data were counted on the following date or dates:

```
all_data %>%
  select(organ, who_plated, who_counted, count_date) %>%
  distinct()
```

```
## # A tibble: 3 x 4
     organ who_plated who_counted count_date
     <chr> <chr>
                       <chr>
                                   <chr>
                                   "\"February 21 2022\""
## 1 lung
            BK
                       BK
## 2 lung
            BK
                       BK
                                   "\"April 18 2022\""
                                   "\"April 25 2022\""
## 3 spleen JR
                       JR
```

#### Distribution of CFUs at each dilution:

#### WE NEED TO ADD SAMPLE CFU PLOTS

Here's a plot that shows how many plates were too numerous to count at each dilution level:

Here is a plot that shows how the CFU counts were distributed by dilution level in the data:

#### 12.5 Identify a good dilution for each sample

# 12.6 Calculate CFUs from best dilution/Estimate bacterial load for each sample based on good dilution

```
# Calculating CFU/ml for every qualifying replicate between 10-75 CFUs. Column binding
meta <- read_excel(path, sheet = "metadata")</pre>
tidy_cfu_meta_joined <- inner_join(tidy_cfu_data, meta) %>%
  group_by(groups) %>%
 mutate(CFUs_per_ml = (CFUs * (dilution_factor^2) * (total_resuspension_mL/volume_pla
 select(organ, count_date, who_plated, who_counted, groups, mouse, dilution, CFUs,
 ungroup()
## Joining, by = "organ"
tidy_cfu_meta_joined
## # A tibble: 146 x 9
     organ count_date
                                who_plated who_counted groups mouse dilution CFUs
##
     <chr> <chr>
                                <chr>
                                           <chr>
                                                      <chr> <chr>
                                                                      <dbl> <dbl>
## 1 lung "\"February 21 2022~ BK
                                                      group~ A
                                                                         3
                                                                               53
## 2 lung "\"February 21 2022~ BK
                                                                         5
                                           BK
                                                      group~ A
## 3 lung "\"February 21 2022~ BK
                                           BK
                                                      group~ A
                                                                         6
## 4 lung "\"February 21 2022~ BK
                                                                          3 119
                                           BK
                                                      group~ B
## 5 lung "\"February 21 2022~ BK
                                           BK
                                                      group~ B
                                                                              48
## 6 lung "\"February 21 2022~ BK
                                           BK
                                                                         5 18
                                                      group~ B
## 7 lung "\"February 21 2022~ BK
                                           BK
                                                      group~ C
                                                                         3 120
## 8 lung "\"February 21 2022~ BK
                                           BK
                                                      group~ C
                                                                              32
## 9 lung "\"February 21 2022~ BK
                                           BK
                                                                          3
                                                                               53
                                                      group~ D
## 10 lung "\"February 21 2022~ BK
                                                                               31
                                                      group~ D
## # ... with 136 more rows, and 1 more variable: CFUs_per_ml <dbl>
```

# 12.7 Create initial report information for these data

#### 12.8 Sample ANOVA

```
cfu_stats <- tidy_cfu_meta_joined %>%
  group_by(organ) %>%
 nest() %>%
  mutate(aov_result = map(data, ~aov(CFUs_per_ml ~ groups, data = .x)),
         tukey result = map(aov result, TukeyHSD),
         tidy_tukey = map(tukey_result, broom::tidy)) %>%
  unnest(tidy_tukey, .drop = TRUE) %>%
  separate(contrast, into = c("contrast1", "contrast2"), sep = "-") %>%
  select(-data, -aov_result, -tukey_result, -term, -null.value)# %>%
## Warning: The `.drop` argument of `unnest()` is deprecated as of tidyr 1.0.0.
## All list-columns are now preserved.
## This warning is displayed once every 8 hours.
## Call `lifecycle::last lifecycle warnings()` to see where this warning was generated.
 # filter(adj.p.value <= 0.05)</pre>
cfu_stats
## # A tibble: 9 x 7
```

```
## # Groups: organ [2]
    organ contrast1 contrast2 estimate conf.low conf.high adj.p.value
    <chr> <chr> <chr> <dbl> <dbl> <dbl>
                                                               <dbl>
## 1 lung group_2 group_1 -15.0 -39.4
                                                   9.34
                                                                0.377
                               -13.1 -39.2 13.1
-2.57 -27.1 22.0
## 2 lung group_3 group_1
                                                                0.562
## 3 lung group_4 group_1
                                                               0.993
## 4 lung group_3 group_2
                                 1.98 -22.7 26.7
                                                              0.997
## 5 lung group_4 group_2
                                12.5
                                         -10.5 35.5
                                                              0.491
                                         -14.4 35.4
## 6 lung group_4 group_3
                                 10.5
                                                               0.689

    -21.5
    -48.8
    5.80
    0.146

    -17.6
    -45.9
    10.7
    0.294

    3.90
    -23.4
    31.2
    0.935

## 7 spleen group_2 group_1 -21.5
## 8 spleen group 3 group 1
## 9 spleen group_3 group_2
```

- 12.9 Save processed data to database
- 12.10 Example one
- 12.11 Example two

# Enzyme-linked immunosorbest assay (ELISA)

ELISA is a standard molecular biology assay for detecting and quantifying a variety of compounds, including peptides, proteins, and antibodies in a sample. The sample could be serum, plasma, or bronchoalveolar lavage fluid (BALF).

#### 13.0.0.1 Importance of ELISA

An antigen-specific reaction in the host results in the production of antibodies, which are proteins found in the blood. In the event of an infectious disease, it aids in the detection of antibodies in the body. ELISA is distinguishable from other antibody-assays in that it produces quantifiable findings and separates non-specific from specific interactions by serial binding to solid surfaces, which is often a polystyrene multiwell plate.

In IMPAc-TB project, it is crucial to evaluate the if the vaccine is eliciting humoral immunity and generating antibodies against vaccine antigen. ELISA will be used to determine the presence of Immunoglobulin (Ig) IgG, IgA, and IgM in the serum different time points post-vaccination.

#### 13.0.0.2 Principle of ELISA

ELISA is based on the principle of antigen-antibody interaction. An antigen must be immobilized on a solid surface and then complexed with an enzymelinked antibody in an ELISA. The conjugated enzyme's activity is evaluated by incubating it with a substrate to yield a quantifiable result, which enables detection. There are four basic steps of ELISA:

- 1. Coating multiwell plate with antigen/antibody: This step depends on what we want to detect the sample. If we need to evaluate the the presence of antibody, the plate will be coated with the antigen, and vice versa. To coat the plate, a fixed concentration of antigen (protein) is added to a 96 well high-binding plate (charged plate). Plate is incubated over night with the antigen at 4 degree celsius (as proteins are temperature sensitive) so that antigens are completely bound to the well.
- 2. Blocking: It is possible that not each and every site of the well is coated with the targeted antigen, and there could be uncovered areas. It is important to block those empty spaces so that primary antibody (which we will add to the next step) binds to these spaces and give us false positive results. For this, microplate well surface-binding sites are blocked with an unrelated protein or other substance. Most common blocking agents are bovine serum albumin, skim milk, and casein. One of the best blocking agents is to use the serum from the organism in which your secondary (detection antibody) is raised. For example, if the secondary antibody is raised in goat, then we can use goat serum as a blocking agent.
- **3. Probing**: Probing is the step where we add sample containing antibodies that we want to detect. This will be the primary antibody. If the antibodies against the antigen (which we have coated) are present in the sample, it will bind to the antigen with high affinity.
- 4. Washing: After the incubation of sample containing primary antibody, the wells are washed so that any unbound antibody is washed away. Washing solution contains phosphate buffer saline + 0.05% tween-20 (a mild detergent). 0.05% tween-20 washes away all the non-specific interactions as those are not strong, but keeps all the specific interaction as those are strong and cannot be detached with mild detergent.
- **5. Detection**: To detect the presence of antibody-antigen complex, a secondary antibody labelled with an enzyme (usually horseradish peroxidase) is added to the wells, incubated and washed.
- **6. Signal Measurement**: Finally to detect "if" and "how much" of the antibody is present, a chromogenic substrate (like 3,3',5,5'-Tetramethylbenzidine) is added to the wells, which can be cleaved the the enzyme that is tagged to the secondary antibody. The color compund is formed after the addition of the substrate, which is directly proportional to the amount of antibody present in the sample. The plate is read on a plate reader, where color is converted to numbers.

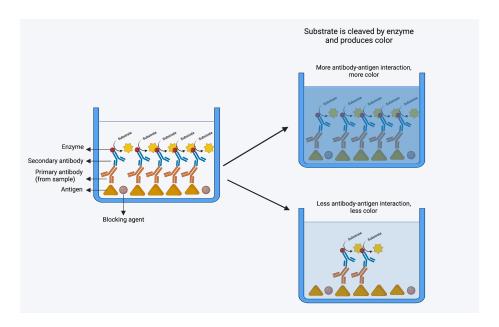


Figure 13.1: A caption

#### 13.1 Read in data from excel file

`` -> `...11` `` -> `...12`

```
library(readxl)
library(tidyverse)

elisa_raw_data <- read_excel("DATA/elisa_s1_07-25-20.xlsx", sheet = "S1", col_names = FALSE, ran

## New names:
## * `` -> `...1`
## * `` -> `...2`
## * `` -> `...3`
## * `` -> `...4`
## * `` -> `...5`
## * `` -> `...6`
## * `` -> `...6`
## * `` -> `...8`
## * `` -> `...8`
## * `` -> `...9`
## * `` -> `...9`
```

#### 13.2 Tidy the data

```
# Convert all columns to numeric

elisa_raw_data_numeric <- elisa_raw_data %>%
    mutate_if(is.character, as.numeric)

## Warning in mask$eval_all_mutate(quo): NAs introduced by coercion

## pivot longer the data

elisa_raw_data_tidy <- pivot_longer(data = elisa_raw_data_numeric, cols = "...1":"...1:

# remove "..." from the first column

elisa_raw_data_tidy$well_id <- str_replace(elisa_raw_data_tidy$well_id, "...", "")

# Add new column to the data_frame

elisa_raw_data_tidy_new <- elisa_raw_data_tidy %>%
    mutate(name = rep(LETTERS[1:8], each = 12))

elisa_raw_data_tidy_new <- elisa_raw_data_tidy_new %>%
```

```
mutate(well_id = paste0(name, well_id)) %>%
  select(-name)
head(elisa_raw_data_tidy_new)
## # A tibble: 6 x 2
   well_id od_450nm
##
    <chr>
                <dbl>
## 1 A1
                0.052
## 2 A2
                0.05
## 3 A3
                0.069
## 4 A4
                0.063
## 5 A5
                0.061
## 6 A6
                0.122
```

#### 13.3 Read in second data set

```
elisa_label_data <- read_excel("DATA/elisa_s1_07-25-20.xlsx", sheet = "S1", col_names = FALSE, n
## New names:
## * `` -> `...1`
## * `` -> `...2`
## * `` -> `...3`
## * `` -> `...4`
## * `` -> `...5`
## * `` -> `...6`
## * `` -> `...7`
## * `` -> `...8`
## * `` -> `...9`
## * `` -> `...10`
## * `` -> `...11`
## * `` -> `...12`
head(elisa_label_data)
## # A tibble: 6 x 12
            . . . 1
                                                    ...2
                                                                            \ldots 3 \quad \ldots 4 \quad \ldots 5 \quad \ldots 6 \quad \ldots 7 \quad \ldots 8 \quad \ldots 9 \quad \ldots 10 \ \ldots 11 \ \ldots 12
##
               <chr>
                                                      <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr> <chr< <chr> <chr
## 1 blank
                                                     secon~ naïv~ 1A-1~ 1A-1~ 1A-1~ 1A-2~ 1A-2~ 1A-2~ 1A-3~ 1A-3~ 1A-3~
## 2 1A-4 (1/250 1A-4 ~ 1A-4~ 1B-1~ 1B-1~ 1B-1~ 1B-2~ 1B-2~ 1B-2~ 1B-2~ 1B-3~ 1B-3~ 1B-3~
## 3 1B-4 (1/250 1B-4 ~ 1B-4~ 2A-1~ 2A-1~ 2A-1~ 2A-2~ 2A-2~ 2A-2~ 2A-2~ 2A-3~ 2A-3~ 2A-3~
```

```
## 4 2B-1 (1/250 2B-1 ~ 2B-1~ 2B-2~ 2B-2~ 2B-2~ 2B-3~ 2B-3~ 2B-3~ 2B-4~ 2B-4~ 2B-4~
## 6 3B-1 (1/250 3B-1 ~ 3B-1~ 3B-2~ 3B-2~ 3B-2~ 3B-3~ 3B-3~ 3B-3~ 3B-4~ 3B-4~ 3B-4~
# pivot longer the data
elisa_label_data_tidy <- pivot_longer(data = elisa_label_data, cols = "...1":"...12",
# remove "..." from the first column
elisa_label_data_tidy$well_id <- str_replace(elisa_label_data_tidy$well_id, "...", "")
# Add new column to the data_frame
elisa_label_data_tidy_new <- elisa_label_data_tidy %>%
 mutate(name = rep(LETTERS[1:8], each = 12))
elisa_label_data_tidy_new <- elisa_label_data_tidy_new %>%
 mutate(well_id = paste0(name, well_id)) %>%
 select(-name)
head(elisa_label_data_tidy_new)
## # A tibble: 6 x 2
## well_id information
##
    <chr>
           <chr>
## 1 A1
           blank
## 2 A2
           secondary
## 3 A3
           naïve (1/250)
## 4 A4
           1A-1 (1/250
## 5 A5
           1A-1 (1/1250
           1A-1 (1/6250
## 6 A6
```

# 13.4 Join the OD table with the information table

```
elisa_data = elisa_raw_data_tidy_new %>% inner_join(elisa_label_data_tidy_new, by="well
head(elisa_data)

## # A tibble: 6 x 3
```

# 13.5 Separate the information tabke into sample ID and dilution columns

```
tidy_elisa_data <- separate(elisa_data, col = "information", into = c("sample_id", "dilution"), s
## Warning: Expected 2 pieces. Missing pieces filled with `NA` in 2 rows [1, 2].
head(tidy_elisa_data)
## # A tibble: 6 x 4
## well_id od_450nm sample_id dilution
## <chr> <dbl> <chr> <chr>
              0.052 "blank"
## 1 A1
                                 <NA>
              0.05 "secondary" <NA>
## 2 A2
             0.069 "naïve " 1/250)
## 3 A3
## 4 A4 0.063 "1A-1 " 1/250
## 5 A5 0.061 "1A-1 " 1/1250
## 6 A6 0.122 "1A-1 " 1/6250
tidy_elisa_data <- tidy_elisa_data %>%
  mutate(dilution = str_extract(dilution, "(/)[0-9]+"),
         dilution = str_replace(dilution, "/", ""),
         dilution = as.numeric(dilution))
tidy_elisa_data <- tidy_elisa_data %>%
  select(well_id, sample_id, dilution, od_450nm)
head(tidy_elisa_data)
## # A tibble: 6 x 4
## well_id sample_id dilution od_450nm
```

#### $40\,CHAPTER\,\,13.\ \ \, ENZYME\text{-}LINKED\,\,IMMUNOSORBEST\,\,ASSAY\,\,(ELISA)$

##		<chr></chr>	<chr></chr>	<dbl></dbl>	<dbl></dbl>
##	1	A1	"blank"	NA	0.052
##	2	A2	"secondary"	NA	0.05
##	3	A3	"naïve "	250	0.069
##	4	A4	"1A-1 "	250	0.063
##	5	A5	"1A-1 "	1250	0.061
##	6	A6	"1A-1 "	6250	0.122

## 13.6 ELISA data analysis optimization

# **Bibliography**

Baazim, H., Antonio-Herrera, L., and Bergthaler, A. (2022). The interplay of immunology and cachexia in infection and cancer. *Nature Reviews Immunology*, 22(5):309–321.